



IsoMet™ High Speed Service Plans



Protect Your Investment

Labs run more efficiently and consistently with well-maintained equipment.




Reduce Unplanned Downtime

Regular equipment maintenance will reduce the hassle of unexpected downtime and associated costs.



Ensure Consistent Results

Confidently analyze lab results with the equipment maintained by the trusted experts at Buehler.

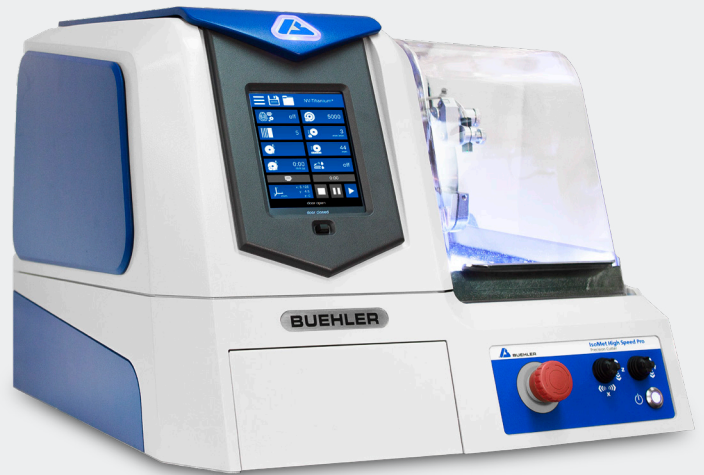

Service Plan Options

Identify the best Service Plan based on the lab's needs.

	Preventative Maintenance	Extended Warranty
Visits Per Year	1	2
Wear Parts Kit	✓	✓
Unlimited Phone Support	✓	✓
Replacement Part Discount	10%	25%
Unlimited Labor Repairs		✓
Priority Response		✓

- Preventative Maintenance ensures consistent results
- Visit includes Wear Parts Kit and 18 Point Maintenance Check list *(see reverse for more details)*
- Extended Warranty must be purchased prior to the expiration of active warranty



Maintenance Checklist

Our comprehensive 18 point maintenance checklist is thorough, ensuring the IsoMet High Speed is running efficiently.

18 Point Maintenance Checklist

- Inspect spindle shaft for alignment
- Check and clean filters on fan
- Inspect interlock switch and hood lock key for proper operation
- Check and clean coolant tank and lines for abrasive builduP
- Check flow of water
- Check vises for proper operation
- Check and inspect coolant lines
- Check front panel components for proper operation
- Check and clean filters on fan
- Check and inspect coolant tank electrical connections
- Inspect drain boot position and condition
- Check and update firmware and software
- Review error log for any potential critical issues with units
- Perform and inspect cuts to ensure machine is working properly
- Inspect hood integrity
- Clean inside cutting area
- Replace hood gas spring




Wear Parts Kit

The parts listed below are included with every Preventative Maintenance visit

IsoMet High Speed Wear Parts Kit 11270SRV

Includes

- Dressing Wheel
- Coolant Tank O-Ring
- Spindle Lock
- Linear Bearing
- Hood Gas Spring
- Coolant Hose

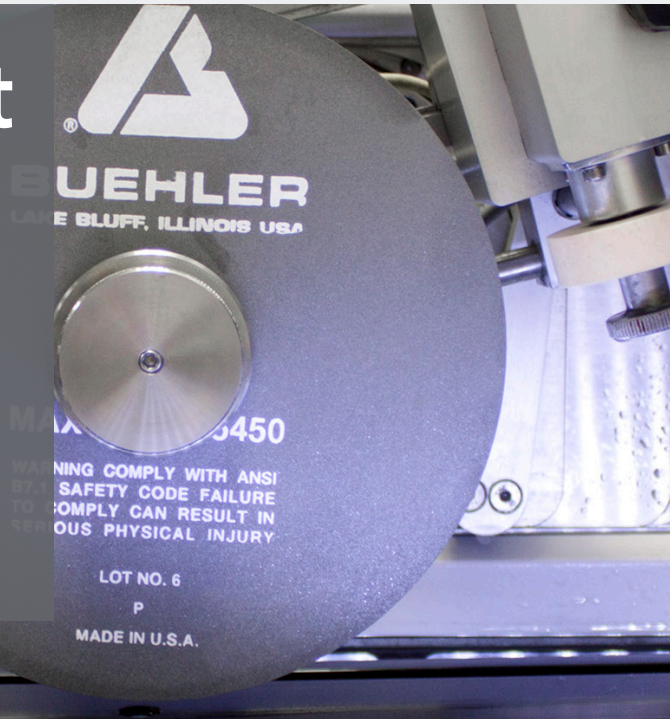

BUEHLER Worldwide Headquarters
North America-South America Offices
41 Waukegan Road
Lake Bluff, Illinois 60044-1699 USA
P: 800 BUEHLER (800-283-4537)
P: (847) 295-6500
www.buehler.com | info@buehler.com

European Headquarters
BUEHLER Germany
info.eu@buehler.com

BUEHLER France
info.eu@buehler.com

BUEHLER United Kingdom
info.eu@buehler.com

BUEHLER China
info.cn@buehler.com

BUEHLER Japan
info.japan@buehler.com

BUEHLER Asia-Pacific
info.asia@buehler.com